POINT-OF-CARE

aptitude Metrix® | Provided by SEKISUI DIAGNOSTICS

Metrix® COVID/Flu Test

A MOLECULAR LAB ANYWHERE AND EVERYWHERE

Features and Benefits

- **Rapid Molecular Test** - Provides lab quality results in 20 minutes
- **3-in-1 Multiplex** - Detects COVID-19, Flu A and Flu B from one swab, in one test
- **Minimal Hands-On Time** - Can be performed in a CLIA-waived setting in a few simple steps and stored at room temperature
- **Convenient** - No maintenance or calibration, ready for testing any time
- **Accessible and Flexible** - Suitable for any facility
- **Exceptional Support** - Experienced support teams and online training modules to help streamline your implementation


Performance Characteristics

Analyte	Positive Percent Agreement*(PPA)	Negative Percent Agreement* (NPA)	Limit of Detection** (LOD)	Time to Result
SARS-CoV-2	96%	99%	500 GE/Swab	20 Minutes
Flu A	96%	99%	1000 GE/Swab	
Flu B	95%	99%	2000 GE/Swab	

*Rounded numbers compared to RT-PCR. See IFU. **Refer to IFU for additional LOD information.





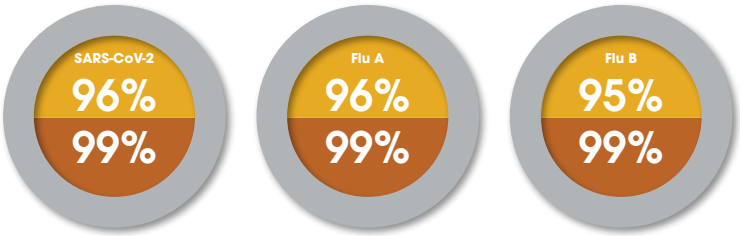

This product has not been FDA cleared or approved, but has been authorized for emergency use by FDA under an EUA; This product has been authorized only for the detection and differentiation of of nucleic acid from SARS-CoV-2, Influenza A, and Influenza B, not for any other viruses or pathogens. The emergency use of this product is only authorized for the duration of the declaration that circumstances exist justifying the authorization of emergency use of in vitro diagnostics for detection and/or diagnosis of COVID-19, Flu A, and Flu B under Section 564(b) (1) of the Federal Food, Drug, and Cosmetic Act, 21 U.S.C. § 360bbb 3(b)(1), unless the declaration is terminated or authorization is revoked sooner.

Metrix® COVID/Flu Test

Experience the Simplicity of Metrix



1. SAMPLE

Break swab in collector
Attach buffer cap
Shake for 20 seconds



2. SNAP

Snap collector onto sensor



3. SLIDE

Slide sensor into reader
Read result within 20 minutes

Interpretation of Results

Results are color-coded for easy interpretation by operators, including co-infection statues[1]



COVID-19 Positive
SARS-CoV-2 viral RNA was detected.



Flu B Positive
Influenza B viral RNA was detected.



Invalid
Test should be repeated.



Flu A Positive
Influenza A viral RNA was detected.



Negative
SARS-CoV-2, Flu A, and Flu B viral DNA not detected.

[1]When more than one virus is detected, the "+" light will flash between the colors of each detected virus

Ordering Information

Test Description	Configuration	Catalog Number	CPT Code
Metrix Reader (Gen 2)	1 Reader	MTRX-RDR-Gen2	
Metrix COVID/Flu Test	25 Tests	MTRX-C19+Flu-25PK	87636QW
Metrix COVID/Flu Molecular Test Bundle	1 Reader + 50 Tests	MTRX-C19+FLU-B1	

To learn more about all of our Metrix Tests, scan the code or go to metrixmolecular.com

Provided by
SEKISUI Diagnostics, LLC
One Wall Street
Burlington, MA 01803
Phone: 800 332 1042
Email: questions@sekisui-dx.com
sekisuidiagnostics.com

Designed and manufactured by
Aptitude Medical Systems Inc.
125 Cremona Drive, Suite 100
Goleta, CA 93117

aptitude | SEKISUI DIAGNOSTICS
Because every result matters™

80-9037-00-00 04/25

©2024 SEKISUI Diagnostics, LLC. All rights reserved. Because every result matters™ is a trademark of SEKISUI Diagnostics, LLC. Metrix® is a registered trademark of Aptitude Medical Systems Inc.