QuickVue™

Influenza+SARS Test

Quick, reliable results for flu and COVID-19 test in only 10 minutes

The QUICKVUE Influenza + SARS Tests can differentiate between SARS-CoV-2, influenza A and influenza B antigens directly from anterior nares swab specimen.


- Simple workflow
- Two-color visual read result
- Actionable results available in as few as 10 minutes

Product performance

Outstanding clinical performance with over 99% of negative percent agreement (NPA) and over 85% of positive percent agreement (PPA).

	SARS-CoV-2	Influenza A	Influenza B
Positive Percent Agreement (PPA)	87.8%	87.2%	87.9%
Negative Percent Agreement (NPA)	99.8%	99.7%	99.7%


Nasal swab procedure

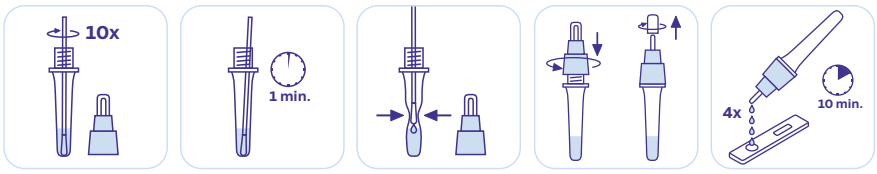

Roll | Leave sample | Squeeze | Add cap | Add sample

Results


Negative **Positive**

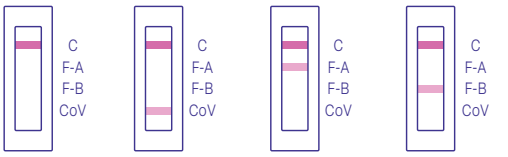

COVID-19 Positive | Influenza A Positive | Influenza B Positive

If the control line at "C" is visible and any other line or multiple lines at 'F-A', 'F-B" and/or 'CoV' appear, the test is positive for that virus.

Invalid

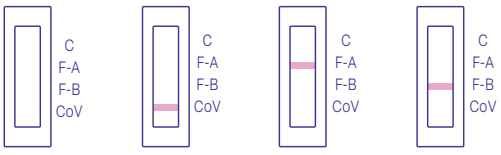

If a line is not visible at 'C', even if any other line is visible in the results window, the result is considered invalid.

Ordering information

Catalog no.	Description	Kit size
20518	QUICKVUE Influenza + SARS Test Kit	25 test kit
20519	QUICKVUE Influenza + SARS Test Controls Kit	

QuidelOrtho™

QuidelOrtho and QUICKVUE are trademarks of QuidelOrtho Corporation or its affiliates and are registered in the U.S. and other countries.
New QuidelOrtho branding may not be available in all markets, subject to country-specific regulatory approval. Please confirm with your local commercial team.
© 2025 QuidelOrtho Corporation | quidelortho.com | 800-874-1517 | PR-109068-CRP-EN-US-v1 (9/25)