Humane Restraint®


PRODUCT SPECS

POLYURETHANE NON-LOCKING BED RESTRAINTS

IMPORTANT PRODUCT INFORMATION—PLEASE READ BEFORE USING



Wrist and ankle bed restraints shown with torso restraint.

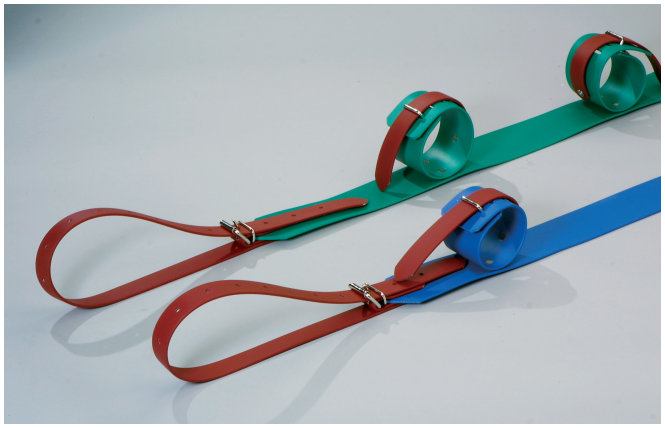

INTENDED USE: To restrict motion of the limbs of the restrained person. Use only as specified in the instructions.







INSTRUCTIONS: Restraints are color coded for quick visual identification. Wristlets are blue and anklets are green. Remember, a person raises their hands to the blue sky and walks on the green grass. Another way to tell is that the wristlets are spaced 26" apart and the anklets are 14" apart.

Lay the restraint flat on the bed with the straps hanging off the side. Loop the straps around the anchoring point and pull through the roller buckle and keeper. Use the same pro-cedure for the strap on the other side.

Then open the restraint by pulling out the polyurethane strap on the wrist or ankle cuff so it is no longer through the loop. Pull up on the end to release from the loop and lay flat. Put the wrist or ankle on the flattened restraint, wrap the cuff around the limb, and adjust to the proper size. Feed the end of the attached strap through the loop to secure.

NOTE: Do not use for prolonged periods of time. Improper adjustment may prohibit proper circulation.

	INVENTORY ID	APPX. WEIGHT
Wrist	CBW-201-POLY	2 lb 2 oz
Ankle	CBA-201-POLY	2 lb 3 oz

DESCRIPTION: The center section is 38" long and 3" wide. The polyurethane side straps are 32" long and 1¼" wide with twelve size adjustments. Wrist cuffs are 12" long and 2¾" wide with four size adjustments ranging from 8" – 10" and ankle cuffs are 14" long and 2¾" wide with four size adjustments ranging from 9½" – 12". Wristlets are spaced at 26" and anklets are 14" apart. Custom sizing available upon request.

SAFETY CHECK: Including but not limited to, loops, rivets, roller buckles, and polyurethane.

MATERIALS USED: Polyurethane, steel loops, and 1¼" steel roller buckles.

CLEANING AND MAINTAINING: For regular cleaning, wipe the polyurethane with a mild soap and warm water. When decontaminating bloodborne pathogens, use one cup of bleach to a gallon of water, then rinse and dry.

STORAGE: Store away from heat greater than 280°F or open flame.

Caution: Federal law restricts this device to sale by or on the order of a physician.

Rev 2/4/05

www.humanerestraint.com ▪ 800.356.7472 ▪ FAX 608.849.6315 ▪ 912 Bethel Circle, Waunakee, WI 53597