
Humane Restraint®

PRODUCT SPECS

LEATHER NON-LOCKING RESTRAINTS


IMPORTANT PRODUCT INFORMATION—PLEASE READ BEFORE USING

INTENDED USE: *To restrict motion of the limbs of the restrained person. Use only as specified in the instructions.*

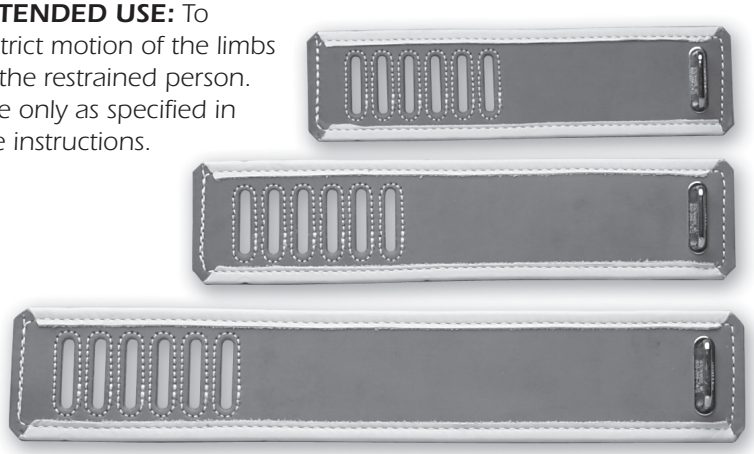




Shown with belt (sold separately).

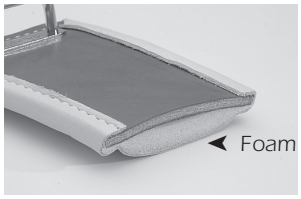

Cross Section

Available with or without inner foam padding for extra comfort.

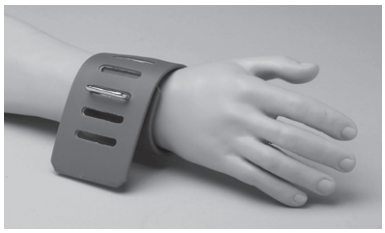



INSTRUCTIONS: *Open restraint completely and lay flat. Put the limb which is to be restrained on the flattened restraint, roll the loop end of the restraint around the limb, then roll the adjustment end over the top of the loop end. Adjust to the proper fitting, allow the loop to go through the appropriate size adjustment, and secure by pulling the tail end of a Humane Restraint belt (locking or non-locking) through the loop and around a stationary object (belt sold separately).*
NOTE: Do not use for prolonged periods of time. Improper adjustment may prohibit proper circulation.

Adjustments are inside circumference in inches and are in increments of ¾".

	WRIST		ANKLE		X-LARGE ANKLE	
	Inventory ID	Approx. Weight	Inventory ID	Approx. Weight	Inventory ID	Approx. Weight
Size Adjustments:	(Six) 5¾"–9½"		(Six) 8"–11¾"		(Six) 11"–14¾"	
Standard	WJ-201	9 oz	AJ -201	10 oz	AJ-201XL	11 oz
Foam Padding	WJL-201	10 oz	AJL-201	11 oz	AJL-201XL	12 oz

DESCRIPTION: *Wristlets are 12" long and 3" wide, anklets are 14" long and 3" wide, and extra large anklets are 17" long and 3" wide.*

SAFETY CHECK: *Including, but not limited to, stitching, loops, and leather.*

MATERIALS USED: *12/13 oz leather, light orthopedic leather lining, steel loops, and ⅜" poly-foam (foam-lined restraints only).*

CLEANING AND MAINTAINING: *For regular cleaning, wipe with a mild soap and warm water. When decontaminating bloodborne pathogens, use one cup of bleach to a gallon of water, then rinse. After wiping the leather, allow it to air dry and apply neat's-foot oil.*

STORAGE: *Store in a dry area to prevent mildew.*

Caution: Federal law restricts this device to sale by or on the order of a physician.

Rev 4/14/05

www.humanerestraint.com ▪ 800.356.7472 ▪ FAX 608.849.6315 ▪ 912 Bethel Circle, Waunakee, WI 53597